

Product Catalogue Number FP045 For Research Use Only – RUO

1p/19q Gene Deletion Probe Detection Kit

[Product Name] 1p/19q Gene Deletion Probe Detection Kit (Fluorescence In Situ Hybridization Method)

[Product Introduction]

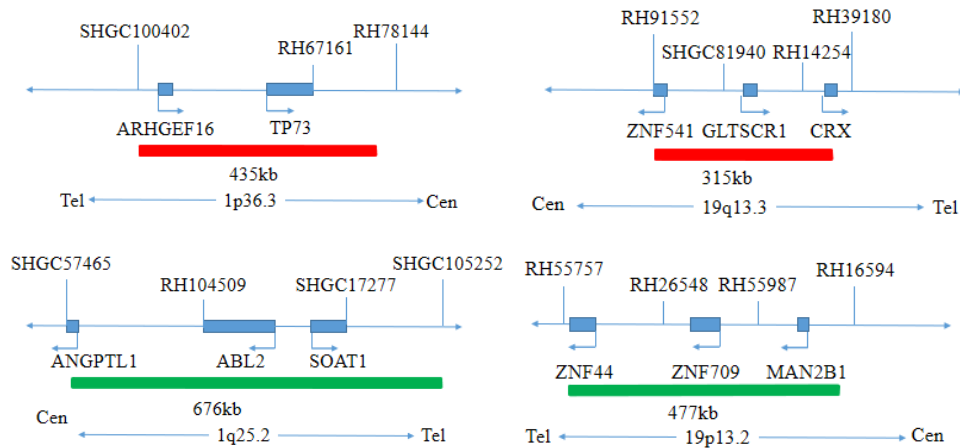
This kit uses orange fluorescent probes 1p36 and 19q13, green fluorescent probes 1q25 and 19p13 to bind 1p/19q probe to the target detection site by in situ hybridization.

[Product Main Components]

The kit consists of 1p36/1q25 dual-color probe and 19q13/19p13 dual-color probe as shown in Table 1.

Table 1 Kit composition

Package	Component name	Specifications	Quantity	Main components
10 Tests/box	1p36/1q25 dual color probe	100µL/Tube	1	1p36 orange probe, 1q25 green probe
	19q13/19p13 dual color probe	100µL/Tube	1	19q13 orange probe; 19p13 green probe



[Storage conditions & Validity]

Keep sealed away from light at -20°C±5°C. The product is valid for 12 months. Avoid unnecessary repeated freezing and thawing that should not exceed 10 times. After opening, within 24 hours for short-term preservation, keep sealed at 2-8°C in dark. For long-term preservation after opening, keep the lid sealed at -20°C±5°C away from light. The kit should be shipped below 0°C.

[Applicable Instruments]

Fluorescence microscopy imaging system including fluorescence microscopy and filter sets suitable for DAPI (367/452), Green (495/517), and Orange (547/565).

[Sample Requirements]

1. Applicable specimen types: Surgical resection or biopsy paraffin-embedded specimens.
2. The tissue should be fixed with 4% neutral formaldehyde fixative solution within one hour after tissue isolation, tissue dehydration and paraffin embedding.

[Instructions]

1. Pre-hybridization treatment

It is recommended to use Wuhan HealthCare Biotechnology Co., Ltd.'s "FISH Pretreatment Reagent" (Catalogue number: CL-003) for pre-treatment operations.

2. Denaturation and Hybridization

The following operations should be performed in a darkroom.

- ① Take the probe at static room temperature for 5 minutes. Briefly centrifuge manually (do not use vortex or shaker instrument). Take 10µl droplet in the cell and drop in the hybridization zone, immediately cover 22mmx22mm glass slide area; spread evenly without bubbles the probe under the glass slide covered area and seal edges with rubber (edge sealing must be thorough to prevent dry film from affecting the test results during hybridization).
- ② Place the glass slide in the hybridization instrument, denature at 85°C for 5 minutes (the hybridizer should be preheated to 85°C) and hybridize at 42°C for 2 to 16 hours.

3. Washing

The following operations should be performed in a darkroom.

- ① Take out the hybridized glass slides, remove the rubber on the coverslip and immediately immerse the slides in a 2xSSC solution for 5 seconds and remove the coverslip.
- ② Place the slides in a 2xSSC at room temperature for 1 min.
- ③ Take out the slides and immerse in a preheated at 68°C 0.3% NP-40/0.4xSSC solution and wash for 2min.
- ④ Remove the slides and immerse in a 37°C preheated deionized water, wash for 1min and dry the slides naturally in the dark.

4. Dyeing

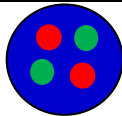
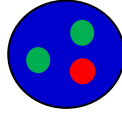
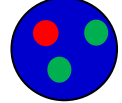
The following operations should be performed in a darkroom

10µl DAPI compound dye is dropped in the hybridization area of the glass slide and immediately covered. The suitable filter is selected for glass slide observation under the fluorescence microscope.

5. FISH results observation

Place the counterstained film under the fluorescence microscope, and first put it under the low-power objective lens (10 ×) Confirm the cell area under the microscope; Go to 40× Under the objective lens, find a position where the cells are evenly distributed; Then in the high-power objective (100 ×) The FISH results of nuclei were observed.

[Common Signal Type Interpretation]

● 1p36 and 19q13 signal	● 1q25 and 19p13 signal
	Negative: 2 Orange-red (2R) ; 2 Green (2G)
	Positive: 1 Orange-red (1R) ; 2 Green (2G) ----1p36/1q25 point-out, 1p36 missing.
	Positive: 1 Orange-red (1R) ; 2 Green (2G) ----- 19p13/19q13 point-out, 19q13 missing.

[Precautions]

- ① This product is for research use only.
- ② The results of this kit will be affected by various factors of the sample itself, but also limited by hybridization temperature and time, operating environment and the limitations of current molecular biology technology, which may lead to wrong results.
- ③ Users must understand the potential errors and accuracy limitations that may exist in the detection process.
- ④ All chemicals are potentially dangerous. Avoid direct contact and wastes should be properly disposed off.

[Manuscript version and approval date]

Manual version: V1.2 reviewed on 07 December 2021

Approval date: 01 April 2019